

Mouse Monoclonal Antibody to

SAPK2 α /p38 α

clone 13D5

Order No.: 0034-100/SAPK2 α -13D5
Size (μ g) 100
Lot No.: 0034S



www.nanotools.de

orders & support:

email: info@nanotools.de

phone: +49-7641-455 670

fax: +49-7641-455 671

05/040613F

Isotype	Species Reactivity	Applications	Mol. Weight	Ref.Cell Line	Epitope	Immunogen
IgG1	human, mouse, rat, dog	WB, ELISA, IP, ICC	38 kDa	SW480	N-terminus, interacts with SAPK2 (p38/CSBP)	peptide conjugated to KLH

Background and Specificity:

Stress-activated Protein Kinases (SAPKs) are strongly activated in response to adverse stimuli such as heat and osmotic shock, UV light and other DNA-damaging reagents, and inhibitors of protein synthesis. They are also activated strongly in response to agonists that are released or produced under conditions of stress, such as proinflammatory cytokines. Stress-activated Protein Kinase 2 (SAPK2/p38) phosphorylates conserved Ser - Pro or Thr - Pro peptide motifs.

Mab SAPK2a-13D5 specifically interacts with the N-terminus of SAPK2 / p38 / CSBP.

Related Products

mab to SAPK1 α /jnk2 (N-terminus)

#0190-100/SAPK1a-12C5

mab to SAPK1 γ /jnk1 (N-terminus)

#0200-100/SAPK1g-5D10

mab to SAPK1/2 (T - G/P - pY)

#0041-100/SAPK1/2-9H8

mab to SAPK2 α (N-terminus)

#0035-100/SAPK2a-20B11

mab to SAPK2 δ (N-terminus)

#0053-100/SAPK2d-5H7

mab to MAPK 1/2 (pT-E-pY)

#0012-100/MAPK-12D4

mab to MAPK 2 (C-terminus)

#0011-100/MAPK2-6G11

mab to MAPK 2 (N-terminus)

#0178-100/MAPK2-6H3

mab to MAPK 2 (internal sequence)

#0239-100/MAPK2-12A4


mab to MAPK7/erk5

#0223-100/MAPK7/erk5-12F2

mab to Mxi 2 (N-terminus)

0046-100/Mxi-2F2

Purification:	The antibody was purified from serum-free cell culture supernatant by subsequent thiophilic adsorption and size exclusion chromatography.
Formulation:	lyophilized from 1 ml PBS / 0.09 % Na-azide / PEG and Sucrose
Reconstitution:	Reconstitute with 1 ml H ₂ O (15 min, RT).
Stability:	For long-term storage, freeze lyophilizate upon arrival (-20°C). Upon reconstitution, aliquote and freeze in liquid nitrogen; reconstituted antibody can be stored frozen at -80°C up to 1 year. Thaw aliquots at 37°C. Thawed aliquots may be stored at 4°C up to 3 months.
	Avoid repeated freeze / thaw cycles.
Positive Control:	#0801: Cell lysate from untreated SW480 cells
Immunoblotting:	1 μ g/ml for HRPO/ECL detection Recommended blocking buffer: Casein/Tween 20 based blocking and blot incubation buffer, e.g. nanoTools product #3031-500/CPPT or #3031-3000/CPPT.
Immunoprecipitation:	use at 1 - 10 μ g per 10 ⁶ pervanadate-treated A431 cells
Immunocytochemistry:	ND
ELISA:	use at 0.05 μ g/ml

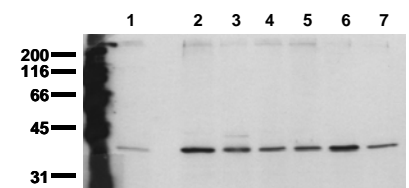

Detection of endogenous SAPK2alpha

Whole cell lysates of serum starved tumor cells (20.000 cells per lane) were applied to SDS-PAGE and transferred to a PVDF membrane. The immunoblot was probed with mab SAPK2a-13D5 (0.5 μ g/ml) for 1h at RT and developed by ECL (exp. time: 30 sec).

lane 1: A431; lane 2: SW480; lane 3: SW620; lane 4: HT29; lane 5: MCF7; lane 6: MDA-MB231; lane 7: T47D

All products are supplied for research and investigational use only. Not for use in humans or laboratory animals.